



**Powerology**

**Therapeutic Massager**  
with Extension Handle Accessory

SKU: PWTMEHA

# Table of Contents

Product Overview.....	2
Features.....	2
Product Heads.....	3
Button Function.....	4
Adjusting the Speed.....	4
Display Interface.....	4
Device Charging.....	4
Specifications.....	5
Warranty.....	6
Contact Us.....	6

Before installing and using the product, please carefully read this User Manual to guarantee correct usage and keep it secure for future reference.

## ◀ Product Overview

Powerology Massage Gun with Extension Handle is a lightweight solution to remove all your muscle tension in no time. It comes with four different heads which make it perfect to use for all parts of your body. The product's base enables wireless charging. In order to fit in all sort of needs it is designed with four different modes as pain relief, calm, energy, and rest, generated with different motor speed from 1500RPM to 3100RPM.

This Product is elegantly crafted for ease of use. This design integrates smoothly into your living space.



It comes with an additional Extension Handle accessory for those hard-to-reach areas.




## ◀ Features


1. 2200mAh Built-in Battery
2. 4 Hours Charging Time
3. 4 Hours Working Time
4. Wireless Charging Base
5. Four Speeds: 3100-1500RPM
6. Less than 55dB Noise Level
7. 1.54" Touch Display

## Product Heads

This product comes with a set of four unique heads, including the dome, cone, thumb, and double. Each is designed to provide targeted therapy for various parts of the body, ensuring comprehensive usage.

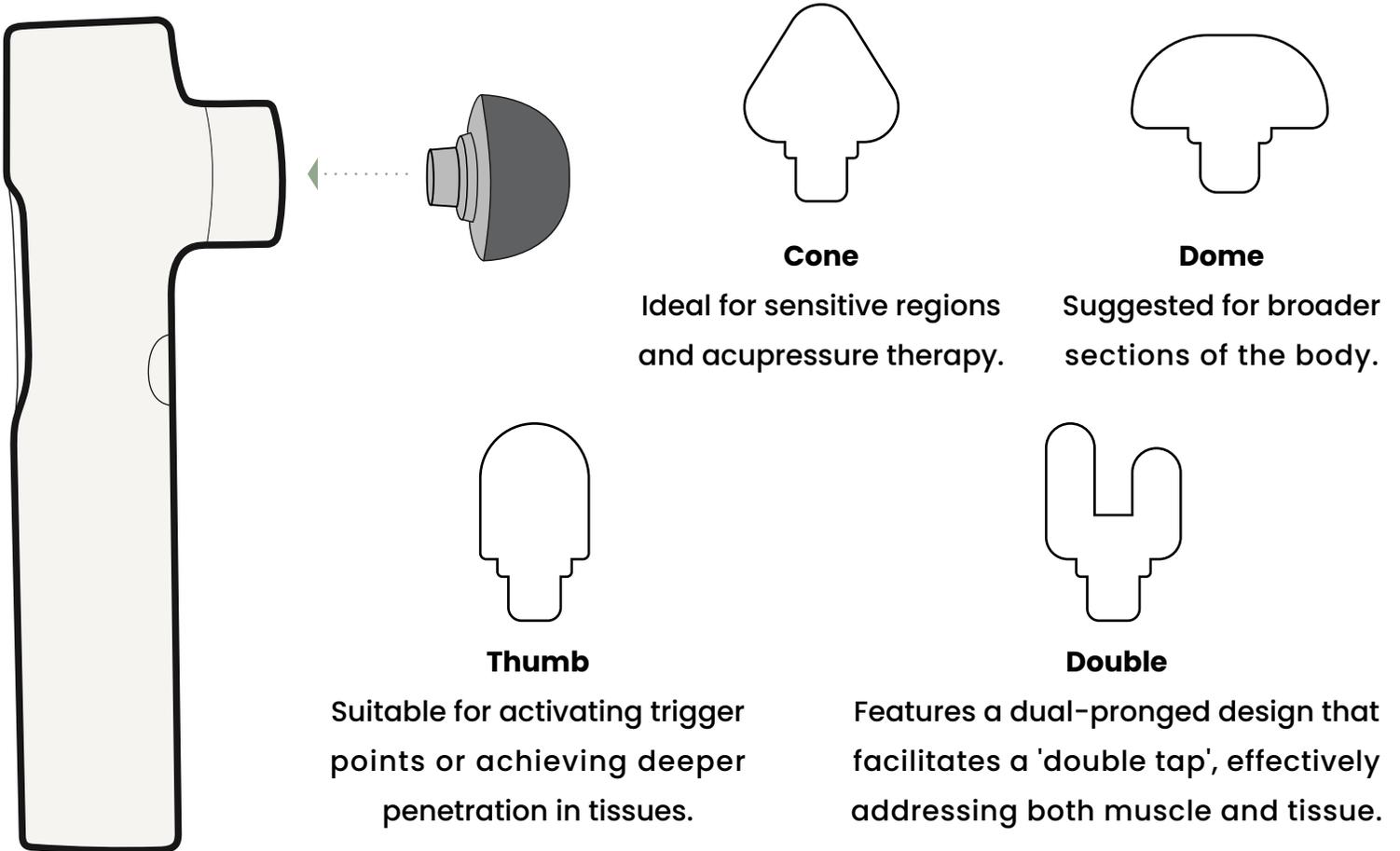

After using the product, the massage head can be conveniently placed back in the base for future use and it facilitates wireless charging for the fascia gun.

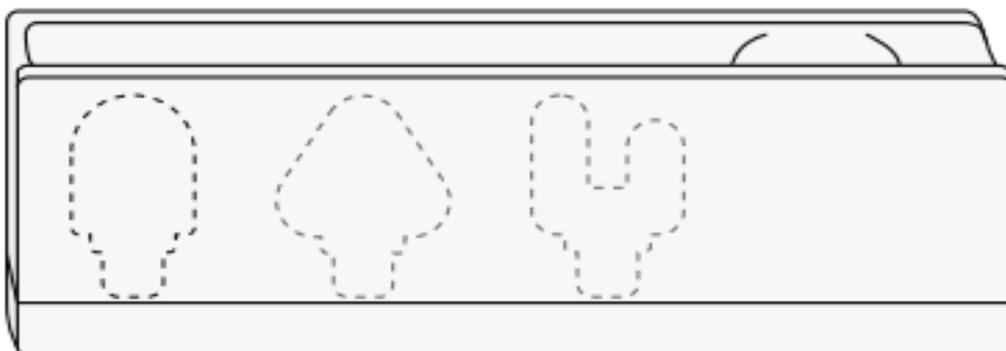


## ◀ Button Function

This product features a single button for powering the device ON or OFF. To activate the device, press and hold the button for 2 seconds. To navigate to the menu, briefly press the button. To go back to the previous menu, press and hold for 4 seconds. For shutting down or resetting the device, press and hold for 8 seconds.


## ◀ Adjusting the Speed

To accommodate a variety of needs, this device offers four distinct modes: "**Pain Relief**", "**Calm**", "**Energy**", and "**Rest**". To navigate through the various vibration speeds, continually press the trigger button until you discover the setting that delivers the utmost comfort for you.

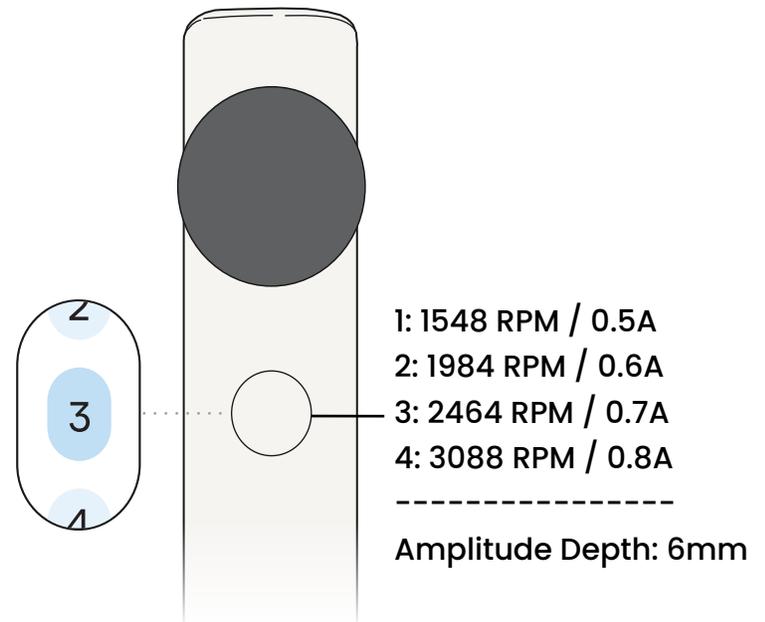

## ◀ Display Interface

Interact with the device's screen by simply swiping your finger across to choose the desired features.

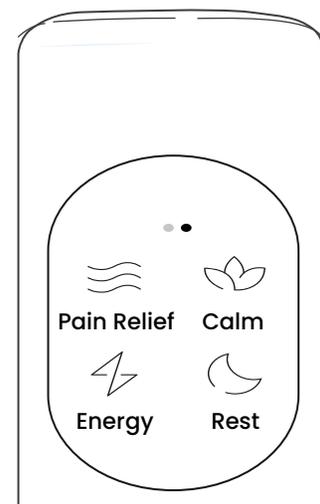

## ◀ Device Charging

This product can be recharged by either of the following methods:

**Method 1:** Connect a charging cable to the Type-C port.

**Method 2:** The docking station is compatible with wireless charging devices.

The charging duration using a Type-C charger is approximately 4 hours, whereas the device can be fully charged using a wireless charger within 16 hours.

**Note:** If the battery level drop below 10 % , the product will engage automatic shutdown protection. The screen may remain active, but the massage function will be inaccessible and the device will power off automatically.

## Specifications

Model	PWTMEHA
Battery Capacity	2200mAh / 7.4V
Battery Type	18650 IP25 Lithium-ion
Power Adapter Input	100-240V, 50/60 Hz
Power	10W
Power Supply	DC 5V/2A
Speed	1500RPM to 3100RPM
Noise	Less than 55dB
Touch Screen Size	1.54"
Motor	Brushless Motor
Charging Port	Type-C port
Discharge Time	≥ 90min
Working Time	About 4 hours (1-2 times a day per week)
Material	ABS+PP+Rubber
Host Size	220×45×72mm
Station Size	222.5×53.5×70mm
Handle Size	220×44×55mm
Host Weight	0.495kg
Station Weight	0.287kg
Handle Weight	0.082kg
Net. Weight	0.864kg
Product Certification	CE

## **Warranty**

Products that you buy directly from our **Powerology** website or shop come with a 24-month warranty.

The 24-month warranty applies to products purchased directly from our **Powerology** website or store. If **Powerology** products are bought from any of our verified retailers, then the product is eligible for only a 12-month warranty. To extend your product's warranty, visit our website **[powerology.me/warranty](https://powerology.me/warranty)** and fill in your details in the provided form along with an uploaded picture of the product to process your request. Once approved, you will receive a confirmation email of the extended product warranty. Upload the required information within 48 hours of purchase to be eligible for a 24-month warranty period.

For more info, please check:  
**[powerology.me/warranty](https://powerology.me/warranty)**

## **Contact Us**

If you have any questions about this Privacy Policy, please contact us at: **[hey@powerology.me](mailto:hey@powerology.me)**

Website: **[powerology.me](https://powerology.me)**

Instagram: **[powerology\\_official](https://www.instagram.com/powerology_official)**

Facebook: **[powerology.ME](https://www.facebook.com/powerology.ME)**